# SANITAS SHK 18

<b>D</b>	<b>Heizkissen</b>	
	<i>Gebrauchsanleitung</i> .....	2
<b>GB</b>	<b>Heating pad</b>	
	<i>Instruction for Use</i> .....	8
<b>F</b>	<b>Coussin chauffant</b>	
	<i>Mode d'emploi</i> .....	14

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**Contents**

1. Items included in the package .....9	4.4 Setting the temperature.....11
1.1 Description .....9	4.5 Automatic switch-off .....11
2. Important instructions	4.6 Switching off.....11
Retain for future use .....9	5. Cleaning and maintenance.....12
3. Intended use.....11	6. Storage.....12
4. Operation.....11	7. Disposal.....12
4.1 Safety .....11	8. What if there are problems? .....13
4.2 Initial use.....11	9. Technical data .....13
4.3 Switching on.....11	

**Explanation of symbols**

	Read the instructions		Machine-washable at 40°C
			Do not wash
			Do not bleach
	Do not insert pins		Do not iron
			Iron at medium temperature
	Do not use folded or rucked		Do not dry-clean
	Not to be used by very young children (0 - 3 years).		The textiles used for this device meet the stringent human ecological requirements of Oeko Tex Standard 100, as verified by Hohenstein Research Institute.
	<b>WARNING:</b> Warning of risks of injury or health hazards		


	<p><b>CAUTION:</b> Safety information about possible damage to appliance/accessories.</p>
	<p><b>NOTE:</b> Important information.</p>

## 1. Items included in the package

- 1 Heating pad
- 1 Cover
- 1 Controller
- 1 Instructions for use

### 1.1 Description

- 1. Mains plug
- 2. Mains cable
- 3. Controller
- 4. Illuminated display for temperature settings
- 5. Sliding switch for ON/OFF and temperature settings
- 6. Plug-in coupling
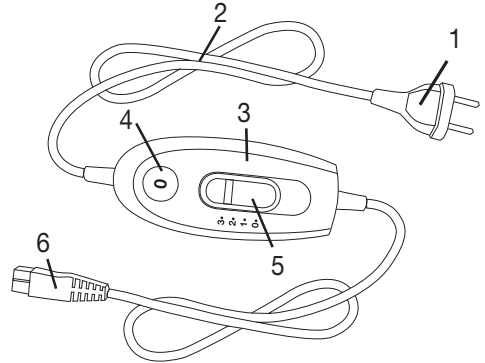
## 2. Important instructions

### Retain for future use

### WARNING

- Non-observance of the following notes may result in personal injury or material damage (electric shock, skin burns, fire). The following safety and hazard information is not only intended to protect your health and the health of others, it should also protect the product. For this reason, pay attention to these safety notes and include these instructions when handing over the product to others.
- This heating pad must not be used by people who are not sensitive to heat or by vulnerable people who may not be able to react to overheating.
- This heating pad may not be used by very young children (0-3 years old) as they are unable to respond to overheating.
- This heating pad may not be used by young children (3-8 years old) unless the controller has been set by a parent or responsible adult and the child has been given adequate instructions on how to use the heating pad safely.

- This heating pad may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised and have been instructed on how to use the heating pad safely, and are fully aware of the consequent risks of use.
- Children must not play with the heating pad.
- Cleaning and user maintenance must not be performed by children unless supervised.

This heating pad is not designed for use in hospitals.

- Do not insert needles.
- Do not use when folded or bunched up.
- Do not use if wet.
- This heating pad may only be used in conjunction with the controller specified on the label.
- The electrical and magnetic fields emitted by this heating pad may interfere with the function of a pacemaker. However, they are still well within the limits: electrical field strength: max. 5000 V/m, magnetic field strength: max. 80 A/m, magnetic flux density: max. 0.1 millitesla. Please consult your doctor and the manufacturer of your pacemaker before using this heating pad.
- Do not pull, twist or make sharp bends in the cables.
- Please check this heating pad frequently for signs of wear and tear or damage.

If any such signs are evident, if the heating pad has been used incorrectly or if it no longer heats up, it must be checked by the manufacturer before being switched on again.

- If the mains connection cable of the heating pad is damaged, it must be replaced by the manufacturer, its customer service department or a similarly qualified person to avoid any risk.
- While the heating pad is switched on, the following conditions must be observed:
  - Objects must not be placed on the pad (e.g. suitcase or washing basket).
  - Heat sources, such as hot water bottles, heating pads or similar, must not be placed on the pad.
- The electronic components in the controller warm up when the heating pad is in use. For this reason, the controller must never be covered or placed on the heating pad when it is in use.
- Always observe the notes on operation (section 4), cleaning and maintenance (section 5) and storage (section 6).
- If you should have any questions about using our devices, please contact Customer Services.

### 3. Intended use

#### IMPORTANT

This heating pad is only designed to warm up the human body.

### 4. Operation

#### 4.1 Safety

#### CAUTION

The heating pad is fitted with a SAFETY SYSTEM. This sensor technology provides protection against overheating across the entire surface of the heating pad with an automatic switch-off in the event of a fault. If the SAFETY SYSTEM has switched off the heating pad, the temperature settings are no longer illuminated when switched on. Please note that for safety reasons, the heating pad can no longer be operated after a fault has occurred and must be sent to the specified service address.


#### 4.2 Initial use

#### IMPORTANT

Ensure that the heating pad will not bunch up or become folded during use.

- First, connect the controller to the heating pad by inserting the plug-in coupling.
- Then plug the mains plug into the socket.
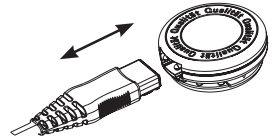
#### 4.3 Switching on

To switch the heating pad on, set the sliding switch for ON/OFF and temperature settings to setting 1, 2 or 3. The temperature settings (4) are illuminated when switched on.

#### 4.4 Setting the temperature

Level 0: OFF

Level 1: minimum heat

Level 2: medium heat

Level 3: maximum heat

#### NOTE

- The fastest way to warm up the heating pad is to initially set the highest temperature setting.

#### WARNING

If the heating pad is being used over several hours, we recommend that you set the lowest temperature setting on the controller to avoid overheating the heated body part, which may lead to burns to the skin.

#### 4.5 Automatic switch-off

This heating pad is equipped with an automatic switch-off function. This turns off the heat supply approx. 90 minutes after the initial use of the heating pad. The temperature setting display then starts to flash. So that the heating pad can be switched back on, the sliding switch for ON/OFF and temperature settings (5) must first be set to setting "0" (OFF). After about 5 seconds it is possible to switch it on again.

## 4.6 Switching off

To switch the heating pad off, set the sliding switch for ON/OFF and temperature settings to setting "0" (OFF). The temperature settings are then no longer illuminated.

### **NOTE:**

If the heating pad is not in use, switch the sliding switch for ON/OFF and temperature settings to setting "0" (OFF) and unplug the mains plug from the socket. Then disconnect the controller from the heating pad by unplugging the plug-in coupling.

## 5. Cleaning and maintenance

### **WARNING**

Before cleaning, always remove the mains plug from the socket first. Then disconnect the controller from the heating pad by unplugging the plug-in coupling. Otherwise there is a risk of electric shock.

### **IMPORTANT**

- The controller must never come into contact with water or other liquids, as this may cause damage.

To clean the controller, use a dry, lint-free cloth. Do not use any chemical or abrasive cleaning agents.

The textile cover can be cleaned in accordance with the symbols on the label and must be removed from the heating pad prior to cleaning.

To clean the plastic protective cover on the heating pad, use a damp cloth.

### **IMPORTANT**

- Please note that the actual heating pad may not be machine-washed, spun or machine dried, chemically cleaned, wrung out, pressed or ironed. Otherwise the heating pad may be damaged.

### **WARNING**

- Only use the heating pad again once the heating pad and cover are completely dry. Otherwise there is a risk of electric shock.

## 6. Storage

If you do not plan to use the heating pad for a long period of time, we recommend that you store it in the original packaging.

### **IMPORTANT**

Please allow the heating pad to cool down before storing. Otherwise the heating pad may be damaged. To avoid sharp folds in the heating pad, do not place any objects on top of it while it is being stored.

## 7. Disposal

Please dispose of the heating pad in accordance with Directive 2002/96/EC on Waste Electrical and Electronic Equipment (WEEE). In case of queries, please contact the municipal authorities responsible for waste disposal in your area.

## 8. What if there are problems?

Problem	Cause	Solution
The temperature settings are not illuminated while - the controller is properly connected to the heating pad - the mains plug is connected to a working socket - the sliding switch is set to 1, 2 or 3	The safety system has switched the heating pad off permanently.	Send the heating pad and the controller for servicing.

## 9. Technical data

See the rating label on the heating pad.

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